

Wave cut scissors are employed in medical settings, particularly for cutting silicone and gel liners used in prosthetics. The serrated, wave-patterned blades help prevent fraying and tearing, ensuring a smoother edge that reduces skin irritation.


Shapes liners for comfortable fit



Wave cut spreads pressure, reducing skin irritation

Usage Instructions:

1. **Clean the Scissors:** Ensure the scissors are sterilized before use to maintain hygiene.
2. **Mark the Material:** Outline the desired cut on the silicone or gel liner.
3. **Steady Cutting:** Hold the material taut and make slow, deliberate cuts along the marked line to maintain the wave pattern's integrity.
4. **Inspect the Edge:** After cutting, check the edge for smoothness and uniformity.

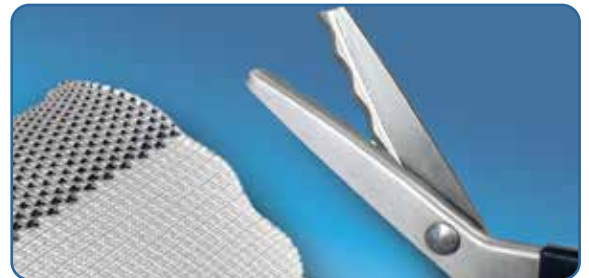

Precautions:

Material Compatibility: Ensure the scissors are suitable for the material you're cutting to avoid damage.

Maintenance: Regularly clean and, if necessary, sharpen the scissors to maintain optimal performance.

Safety: Always handle scissors with care, keeping fingers away from the blades during use.

